
wellion


BONA HbA1c ANALYZER



CLEANING AND MAINTENANCE GUIDE

Overview of analyzer

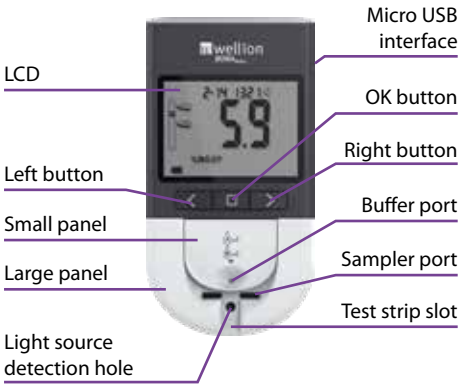

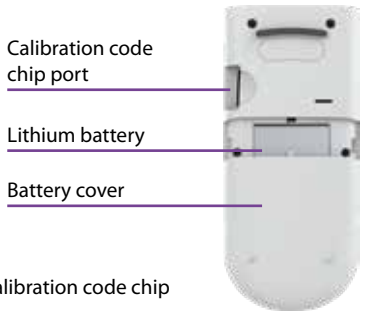


Calibration code chip

Battery installation

1. To open the battery door, slide the battery cover in the direction of the arrow.
2. Place the battery into the battery holder. Make sure the battery is in the correct direction.
3. Put the battery cover back in place.



Cleaning and maintenance

Cleaning operation should be conducted in the following cases.

1. Residues of buffer A, buffer B or blood stains in the analyzer.
2. After daily use.
3. Analyzer contaminated by dust, liquid or any other contaminant.

1.



Analyzer should be turned off before cleaning.

2. Use cotton swab to absorb possible buffer residual in the buffer port.



3. Remove used test strip.




4. Take off the small and large panels together.



5. Clean the surface of panels.



6. Clean the light source area by using cotton swab with distilled water (use alcohol cotton swab if necessary). Please use a dry cotton swab to dry the area.

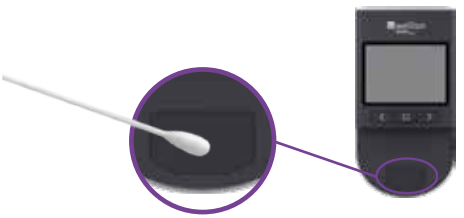


7. Reinstall all panels after cleaning.



If liquid penetrates into analyzer, disassemble the small and large panels and clean the surfaces of both panels.

Disassemble the small and large panels.

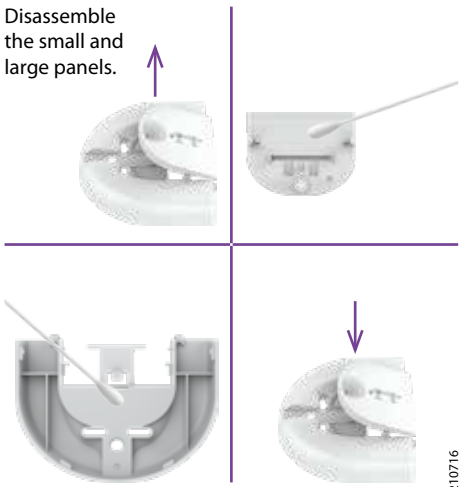

MED TRUST Handelsges.m.b.H.
Gewerbepark 10
7221 Marz, Austria

CE 0197

W16-38EN 20210716